



# Gluteal implants

**ROUND, BICONVEX, ANATOMICAL**



Sebbin  
PARIS

The partner your body can trust



**SEBBIN IS ONE OF THE FIRST EUROPEAN MANUFACTURERS TO OFFER THIS IMPLANTABLE MEDICAL DEVICE.**

We use biocompatible, medical grade, long-term implantable raw materials.

Our manufacturing processes and their monitoring are recorded and stored for a period of 30 years. Our quality system complies with NF EN ISO 9001: 2015 and NF EN ISO 13485: 2016 standards and enables the manufacture of products that comply with European regulations.

Gluteal implants prefilled with silicone gel are more flexible, easier to introduce and offer a better projection compared to solid silicone implants<sup>(1,2,3)</sup>.

A specific very highly cohesive silicone gel was selected so that our gluteal implants rapidly return to their initial shape following the strong stress that occurred at the time of introduction.

Sebbin offers 3 shapes of smooth gluteal implants, to suit the different existing criteria.

**Sources:**

- (1) Some comments on gluteal reshaping. Gonzalez R. 2008.
- (2) Submuscular gluteal augmentation: 17 years of experience with gel and elastomer silicone implants. Hidalgo JE. 2006.
- (3) Subfascial gluteal augmentation. De la Peña JA, Rubio OV, Cano JP, Cedillo MC, Garcés MT. 2006.
- (4) Classification system for gluteal evaluation. Mendieta CG. 2006.




### Round gluteal implants



This shape can be used when:<sup>(4)</sup>

- Most of the gluteal volume is in the lower part of the buttock,
- The gluteal muscle is short,
- The buttocks form a small square.

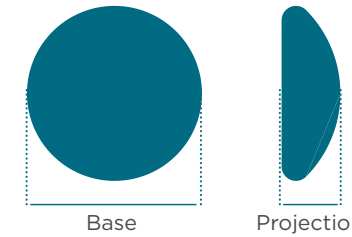


Reference	Volume (ml)	Base (mm)	Projection (mm)
LS 04 200	200	108	30
LS 04 235	235	113	33
LS 04 270	270	120	34
LS 04 300	300	125	36
LS 04 330	330	127	38
LS 04 360	360	128	42
LS 04 400	400	132	43
LS 04 440	440	136	44

### Biconvex gluteal implants



This shape corresponds to the following criteria:

- The biconvex-shaped is for the same indications as the round-shaped but ensures more projection for an identical base.
- In the event of a shift of the implant, the symmetrical-shaped makes this shift imperceptible to the patient.

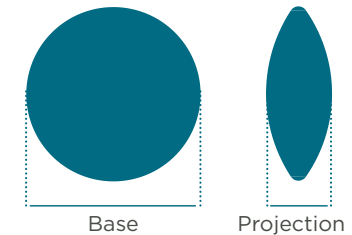

Reference	Volume (ml)	Base (mm)	Projection (mm)
LS 06 370	370	127	52
LS 06 410	410	130	53
LS 06 480	480	136	56
LS 06 530	530	143	59


### Anatomical gluteal implants



This shape can be used when:<sup>(4)</sup>

- Most of the gluteal volume is in the upper part of the buttock,
- The gluteal muscle is long,
- The buttocks form a high square,
- The patient is tall,

The main indication is to increase the volume of the lower part of the buttock.

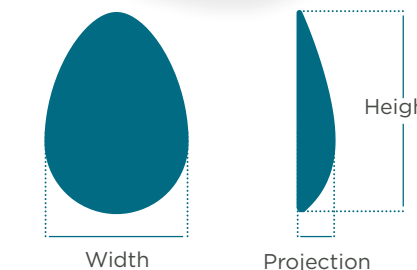

Reference	Volume (ml)	Width (mm)	Height (mm)	Projection (mm)
LS 05 185	185	97	139	30
LS 05 260	260	107	148	32
LS 05 335	335	117	162	36
LS 05 400	400	129	174	42



The partner your body can trust

**Sebbin Benelux SPRL**

Tél.: +32 2 792 57 76 - Fax: +32 2 792 57 73  
Email: info@sebbin.be

**Sebbin Deutschland GmbH**

Tél.: +49 2102 70 88 00-23 - Fax: +49 2102 70 88 00-99  
Email for Germany: service@sebbin.de  
Email for Austria: austria@sebbin.de

**Groupe Sebbin SAS, France**

Tél.: +33 1 34 42 13 28 - Fax: +33 1 34 42 16 88  
Email: serviceclients@sebbin.com

**Sebbin Ibérica Productos Médicos, S.L.**

Tél.: +34 93 159 3003 - Fax: +34 93 676 3053  
Email: info@sebbin.es

**Sebbin Suisse Sarl**

Tél.: +41 43 336 72 30 - Fax: +41 43 336 72 36  
Email: info@sebbin.ch

**Sebbin UK Ltd**

Tél.: 0800 652 4028 - Fax: +4433 0088 9575  
Email: custserv@sebbin.uk

This document is intended for healthcare professionals. According to the Medical Device Directive 93/42/EEC, gluteal implants are class IIb devices, manufactured by Groupe SEBBIN and intended to be used in the context of plastic, reconstructive and aesthetic surgery. These devices are CE marked by the notified body mdc medical device certification GmbH number 0483. Please read the instructions for use carefully before use. They are available on [www.mysebbin.com](http://www.mysebbin.com). These devices cannot be reimbursed by health insurance companies. Our sterile, single-use gluteal implants are sold individually. Dimensions of our round, biconvex and anatomical implants may present slightly variations to the ones mentioned into this document.

[www.sebbin.com](http://www.sebbin.com)

GRUPE SEBBIN SAS

39 Parc d'Activités des Quatre Vents - CS80010 Boissy l'Aillerie - 95651 CERGY-PONTOISE CEDEX  
Tel.: + 33 1 34 42 13 28 - Fax: + 33 1 34 42 16 88